

Beaver[®] Safety Knives for MICS

Innovative Solutions for Micro-Incision Cataract Surgery

Beaver-Visitec is proud to offer a broad range of safety engineered, single-use knives for micro-incision cataract surgery. All Beaver[®] Xstar[®] Safety Knives are designed for exceptional sharpness and consistency.



Beaver [®] Xstar [®] Safety Slit Knife		
Description	Single Bevel	Double Bevel
1.4 mm, angled 45 degrees	378214	-
1.5 mm, angled 45 degrees	378215	-
1.7 mm, angled 45 degrees	378217	-
1.8 mm, angled 45 degrees	378218	-
2.2 mm, angled 45 degrees	378222	378822
2.4 mm, angled 45 degrees	378224	378824
2.5 mm, angled 45 degrees	370145	378825

Beaver [®] EdgeAhead [™] Safety Sideport MVR	
Description	Item #
.90 mm (20G)	378231
.60 mm (23G)	378233

DESIGNED TO MEET
OSHA
SAFETY REQUIREMENTS


For ordering information call:
1.866.906.8080/1.781.906.8080
or visit www.beaver-visitec.com

bvi Beaver Visitec
Keeping Your Vision in Sight

Beaver® Safety Knives for MICS

Safety Engineered to Protect You and Your Staff

Beaver® Safety Knives are designed to minimize the risk of accidental surgical blade injuries that could occur during all phases of a surgical procedure: before use, during passing, between steps of the procedure, after use or during disposal. The protective shield is easily activated by the user through a spring-assisted slider that provides tactile and audible feedback in both the shielded and unshielded positions.

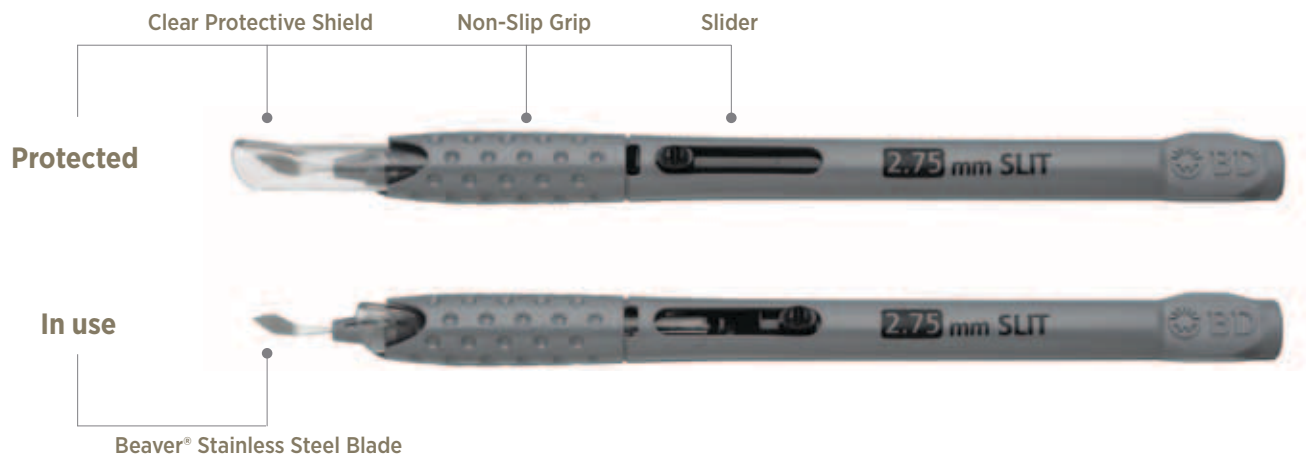

To Use

- Pull slider on handle to retract protective shield.
- Retract shield to rear locked position to fully uncover blade.

To Cover

- Push slider on handle to engage protective shield.
- Advance shield to forward locked position at lock indicator line.

To Dispose of

- Dispose of protected knife in an approved sharps container



Beaver-Visitec International, Inc.

411 Waverley Oaks Road, Waltham, MA 02452

Tel: 1.866.906.8080 Fax: 1.866.906.4304

www.beaver-visitec.com